

DOG MOBILITY AID

Falcon

Dog Wheelchair for Back Legs, Adjustable Frame Rated to 143 lb

Dogs up to 143 lb with hind-leg weakness or paralysis and a strong front end. If the front legs are going too, a rear-support frame is the wrong one.



SPECIFICATIONS


PRODUCT TYPE	Dog Mobility Aid
TYPE	Rear-support, two wheels
FOR	Hind-leg weakness or paralysis, good front-end strength
WEIGHT CAPACITY	143 lb (65 kg)
FRAME	Aircraft-grade aluminum alloy
ADJUSTS	Length, height and width
HARNESS	Padded mesh body sling with chest straps
WHEELS	Spoked, solid tire
REAR PAWS	Stirrup loops, up or down
SURFACES	Grass, gravel, paths and indoors
ASSEMBLY	Frame arrives part-built, fitted to your dog

Dimensions are approximate. Natural materials vary slightly from piece to piece.

KEY FEATURES

- **Adjusts three ways.** Length, height and width, so it is fitted to your dog rather than close enough.
- **Aircraft-grade aluminum.** Rated to 143 lb and light enough that he carries it rather than drags it.
- **Stirrups up or down.** Rear paws can ride or trail, depending on what he can still use.

Dimensions

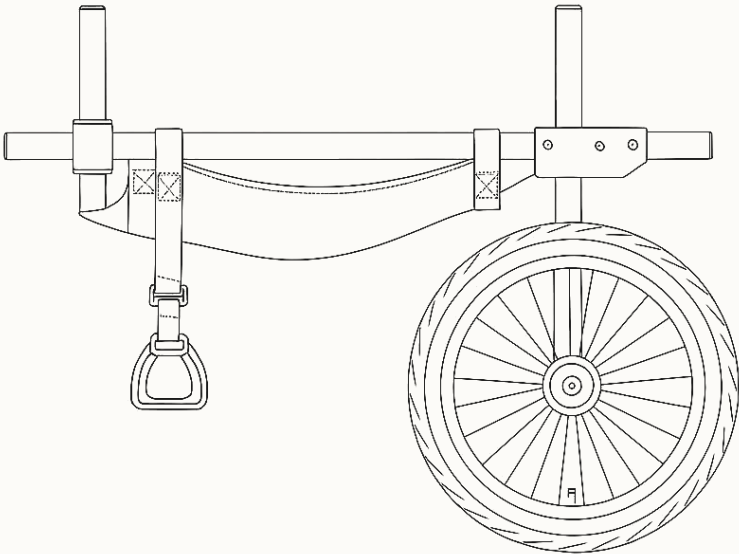

OVERALL VIEW · NOT TO SCALE · DIMENSIONS AS LISTED

MEASUREMENTS

See the specifications on page 1.

VIORN	
ITEM Falcon · Dog Mobility Aid	
DRAWING NO. VI-MA-41730-D1	SCALE Not to scale
UNITS in / cm	REVISION A
ISSUED 9 Oct 2026	

Ordering & details

Catalogue numbers are Viorn reference numbers. Quote the catalogue number when ordering or contacting us.

CONFIGURATIONS

CONFIGURATION	CATALOGUE NO.
Standard	VI-MA-41730-01

WHO IT FITS

Dogs up to 143 lb with hind-leg weakness or paralysis and a strong front end. If the front legs are going too, a rear-support frame is the wrong one.

WILL IT FIT?

You need three measurements before you order: floor to the top of the hips, hip width, and shoulder to the base of the tail. The frame adjusts in all three directions, but it has limits – send us the numbers and a photo and we will tell you whether this chair covers your dog before you buy it. A frame that nearly fits will rub.

MATERIALS

Aircraft-grade aluminum alloy

CARE

Hand-wash or wash as the label directs; air dry. Check straps, buckles and stitching before each use.

WARRANTY & RETURNS

Frames and hardware: 5 years. Covers, foam, fabric and mesh: 1 year. Cover runs from delivery, for the original buyer; chewing, clawing, overloading, misuse and normal wear are not covered. Returns within 30 days of delivery, unused and in original packaging, with no restocking fee. If it arrives damaged or wrong, we pay the return shipping.

Fitting & use

GETTING STARTED

1. **Measure before you order.** Hip height, hip width, shoulder to tail. Send them with a photo. We would rather talk you out of the wrong chair than process a return.
2. **Fit it empty, then fit it on him.** Set the rough adjustment with the frame on the floor, then fine-tune with him standing in it. His spine should stay level - not tipped up at the back.
3. **Ten minutes on grass, not a walk.** Short sessions, you alongside, somewhere soft. Most dogs have it by the third go. Build up from there.

BEFORE YOU START

- ☐ Your dog measured as the size chart directs
- ☐ A vet consulted first for post-surgery or neurological conditions
- ☐ A short, calm first session indoors

PLACEMENT GUIDELINES

FIT

Fit snugly enough not to slip, loose enough to slide two fingers under every strap.

FIRST USES

Start with short sessions and check for rubbing at the armpits, groin and belly after each one.

LIFTING

Lift with your legs, not your back, and take your dog's weight gradually.

Safety. Never leave a dog unattended in a mobility aid. Stop and check the fit if you see rubbing, swelling or distress, and follow your vet's advice.

Questions about sizing or setup?

support@viorn.com · +1 (813) 773-3585

9am – 5pm Eastern, seven days. We reply within 12 hours on business days.

Returns. 30 days from delivery, no restocking fee.