

DOG MOBILITY AID

Osprey

Full-Body Dog Lift Harness, Padded Handle, Chest 25–39"

Dogs who need lifting as a whole body – hind-leg weakness, post-operative recovery, degenerative myelopathy. A 25–32 in chest covers most retrievers and collies; 33–39 in suits a big labrador, shepherd or rottweiler.



SPECIFICATIONS

PRODUCT TYPE	Dog Mobility Aid
TYPE	Full-body lift harness
SIZES	1: chest 25 to 32 in. 2: chest 33 to 39 in
LIFT	A padded panel under the chest and belly, up to one padded handle
WEBBING	Reflective
LINING	Breathable mesh
BUCKLE	Side-release
FOR	Hind-leg weakness, post-operative recovery, degenerative myelopathy, dogs who can stand but cannot get started
NOT FOR	A dog with no hind-leg function at all, who needs a wheelchair

Dimensions are approximate. Natural materials vary slightly from piece to piece.


OPTIONS

SIZE
Chest 25-32 in, Chest 33-39 in

KEY FEATURES

- **Lifts the whole dog.** A padded panel under the chest and belly, so he is held level rather than folded.
- **One padded handle.** The weight comes up through your arm and your back stays straight.
- **Breathable mesh inside.** For a dog who wears it several times a day for weeks.

Dimensions

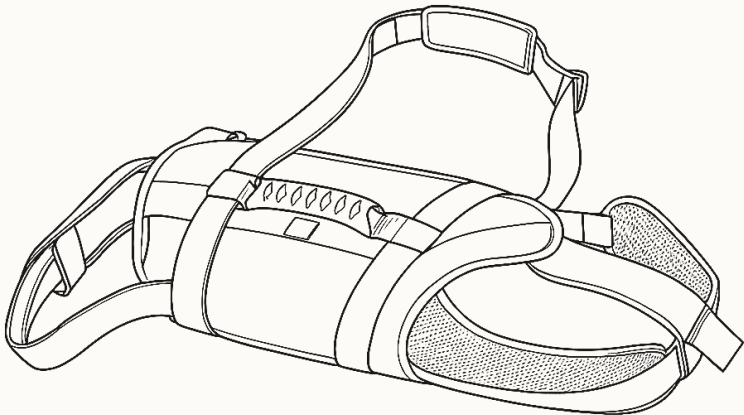

OVERALL VIEW · NOT TO SCALE · DIMENSIONS AS LISTED

MEASUREMENTS

SIZES	1: chest 25 to 32 in. 2: chest 33 to 39 in
-------	--

VIORN	
ITEM Osprey · Dog Mobility Aid	
DRAWING NO. VI-MA-43426-D1	SCALE Not to scale
UNITS in / cm	REVISION A
ISSUED 9 Oct 2026	

Ordering & details

Catalogue numbers are Viorn reference numbers. Quote the catalogue number when ordering or contacting us.

CONFIGURATIONS

SIZE	CATALOGUE NO.
Chest 25-32 in	VI-MA-43426-01
Chest 33-39 in	VI-MA-43426-02

WHO IT FITS

Dogs who need lifting as a whole body — hind-leg weakness, post-operative recovery, degenerative myelopathy. A 25–32 in chest covers most retrievers and collies; 33–39 in suits a big labrador, shepherd or rottweiler.

WILL IT FIT?

Measure the chest at its widest point, just behind the front legs. 25 to 32 inches is most retrievers and collies; 33 to 39 inches suits a big labrador, shepherd or rottweiler. If you are between sizes, take the larger: the straps pull in a long way.

CARE

Hand-wash or wash as the label directs; air dry. Check straps, buckles and stitching before each use.

WARRANTY & RETURNS

Frames and hardware: 5 years. Covers, foam, fabric and mesh: 1 year. Cover runs from delivery, for the original buyer; chewing, clawing, overloading, misuse and normal wear are not covered. Returns within 30 days of delivery, unused and in original packaging, with no restocking fee. If it arrives damaged or wrong, we pay the return shipping.

Fitting & use

GETTING STARTED

1. **Measure the chest first.** At its widest point, just behind the front legs. Take the larger size if you are between two.
2. **Check the skin daily.** Take it off overnight and check the skin under the padding daily, the way you would with any harness used for weight.
3. **Lift with your legs.** You can lift him into the car with it, but keep him close to your body. For a repeated daily lift, a ramp or steps will do your back more good.

BEFORE YOU START

- ☐ Your dog measured as the size chart directs
- ☐ A vet consulted first for post-surgery or neurological conditions
- ☐ A short, calm first session indoors

PLACEMENT GUIDELINES

FIT

Fit snugly enough not to slip, loose enough to slide two fingers under every strap.

FIRST USES

Start with short sessions and check for rubbing at the armpits, groin and belly after each one.

LIFTING

Lift with your legs, not your back, and take your dog's weight gradually.

Safety. Never leave a dog unattended in a mobility aid. Stop and check the fit if you see rubbing, swelling or distress, and follow your vet's advice.

Questions about sizing or setup?

support@viorn.com · +1 (813) 773-3585

9am – 5pm Eastern, seven days. We reply within 12 hours on business days.

Returns. 30 days from delivery, no restocking fee.